

## **Zwick**Materials Testing

### **Product Information**

Vibrophores Amsler 50 – 250 HFP 5100



Figure: Amsler 150 HFP 5100

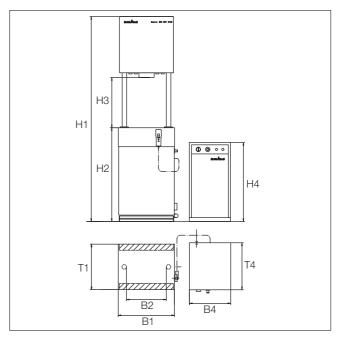


Figure: Drawing of Amsler 10 HFP 5100 with dimensions

#### Range of application

For fatigue testing of materials or components by applying sinusoid loads Zwick Vibrophores offer an economical alternative using the resonance principle with constant or variable amplitude and mean load.

Fatigue testing of materials and components, e.g. fatigue tests according to DIN 50100 (S/N curve) for tensile stress, compressive and alternating stress ranges. Testing can be performed either force controlled or strain controlled.

Typical applications are fracture toughness tests on Compact Tensile (CT) and Single Edge Bending (SEB) specimens, oscillating of bending specimens, fatigue strength tests on components such as springs and screws. Other applications include production and quality control of components that are exposed during their lifetime.

Tests can also be carried out under various environmental conditions (temperature, aggressive media) with corresponding special features. In addition torsion and bending tests can also be carried out.

#### **Advantages and features**

- High specimen throughput due to high test frequencies and short test time
- Extremely low energy consumption due to resonance drive (2% of the cost of a typical servo-hydraulic testing machine)
- Simple installation, no additional power packs necessary (e.g. hydraulics, cooling water)
- Virtually maintenance-free system through lack of wearing parts
- Precision force measurement through calibration according to DIN 51 221 and US MIL Std. 1312 B
- High testing precision through adaptive control
- Fast and precise mean load control by servo drive
- High system stability and low noise-sensitivity through pulse width modulator
- Two measurement channels for additional extensometers (force, strain), optional extensible through two measuring inputs
- User friendly operation through the use of industry standard testXpert® application software configured precisely for each test
- Standard PC can be used because no special interface card is necessary
- The testing frequency can easily be changed using standard weights



# **Zwick**Materials Testing

## **Product Information**

Vibrophores Amsler 50 - 250 HFP 5100

Fmax [kN]	50	100	150	200	250
Mean load max. [kN]	± 50	± 100	± 150	± 150	± 150
Force amplitude max. [kN]	± 25	± 50	± 75	± 100	± 125
Oscillation range of specimen max. [mm]	4	6	6	4	4
Testing frequency range [Hz]	35 – 300	35 – 300	35 – 300	35 – 300	35 – 300
Number of frequency stages <sup>1</sup>	5	5	5	5	5
Dimensions load frame					
H1 – Height [mm]	2700	2700	2700	2700	2700
B1 – Width [mm]	750	750	750	750	750
T1 – Depth [mm]	600	600	600	600	600
Dopar [rimi]	200	300	300	300	300
Dimensions test area					
H2 – Height of machine table [mm]	1235	1235	1235	1235	1235
H3 – Height between machine table and					
load cell [mm]	190 – 690	180 – 680	180 – 680	160 - 660	160 - 660
B2 – Horizontal daylight [mm]	530	530	530	530	530
Travel of machine frame max. [mm]	500	500	500	500	500
Measurement and control electronics VibroWin					
H4 – Height [mm]	1040	1040	1040	1040	1040
B4 – Width [mm]	550	550	550	550	550
T4 - Depth [mm]	620	620	620	620	620
Cable length between HFP and VibroWin [m]	5	5	5	5	5
J					
Weight approx. [kg]	2900	3000	3000	3000	3000
Item no.					
Standard type	• 940738	• 940739	• 940740	• 940750	• 940751
	(BRA350008011)	(BRA350008012)	(BRA350008013)	(BRA350008014)	(BRA350008015)
With column extension 2 150 mm	• 924508	• 924513	• 924516	• 924518	• 924523
	(BRA350008011.A)	(BRA350008012.A)	(BRA350008013.A)	(BRA350008014.A)	(BRA350008015.A)
With column extension <sup>2</sup> 250 mm	• 924509	• 924514	• 924517	• 924519	• 924524
	(BRA350008011.B)	(BRA350008012.B)	(BRA350008013.B)	(BRA350008014.B)	(BRA350008015.B)
With column extension <sup>2</sup> 350 mm	• 924512	• 924515	• 924510	• 924522	• 924525
	(BRA350008011.C)	(BRA350008012.C)	(BRA350008013.C)	(BRA350008014.C)	(BRA350008015.C)

<b>Environmental conditions</b>	
Operating temperature [°C]	+10 +35
Storage temperature [°C]	-25 +55
Humidity range (not condensing) [%]	< 90
Electrical connection	
Mains voltage 3 Ph/N/PE <sup>3</sup> [V]	400
Mains frequency [Hz]	50 / 60
Drive power [kVA]	1
Fuse [A]	16
Noise level at 1m distance [dB(A)]	80 – 110

<sup>&</sup>lt;sup>1</sup> by activation / deactivation of weights

<sup>&</sup>lt;sup>2</sup> with column extension the stated machine heigt and weight of the vibrophore changes

<sup>&</sup>lt;sup>3</sup> Three phase AC motor (L1, L2, L3), neutral wire N, protective earth PE